



User Guide

Corsano 5G Hub — onboarding with the Corsano IoT Hub desktop
From plugging in the hub to vital signs in the Corsano Health Cloud.

Version: 1.1

Date: 2026-08-29

Applies to: Desktop app v2.0.0 — hub firmware $\geq$ 1.33



Table of Contents

1. Overview
2. Prerequisites
3. Onboarding step by step
4. End-to-end verification
5. Troubleshooting
6. Day-to-day operations
7. Support



1. Overview

1.1 How it works

Corsano hub connects to several 287-2B bracelets in its range, captures the vital parameters and uploads the data to the Corsano Health Cloud over WiFi, Ethernet or 5G.

Bracelets 287-2B → Corsano 5G Hub – WiFi / Ethernet / 5G → Corsano Health Cloud

The desktop app is only used for onboarding: register the hub in the cloud, attach it to a study, configure its network. Once provisioned the hub runs on its own — unplug the Computer cable and connect the hub to a power source.

1.2 The hub at a glance



Figure 1 — The Corsano Hub.

Element	Purpose
LED status ring	The hub's only display — shows boot, configuration and operating states.
WPS button	Short press joins a WiFi network via WPS push-button pairing.
USB-C port	Power, and the configuration link to the desktop console.
RJ45 Ethernet	Wired uplink (preferred when available).
Nano-SIM slot	5G cellular uplink — SIM factory-installed.
microSD slot	Service use by Corsano technicians only.



1.3 Security

- Vitals are encrypted and authenticated on the bracelet itself before they are ever sent. The BLE link carries no readable data.
- Everything stored on the hub is encrypted, and the hub only runs firmware signed by Corsano (secure boot).
- The hub carries a dedicated hardware crypto core — a Microchip ATECC608A secure element (Common Criteria EAL4+ certified).
- All uplinks to the Corsano Health Cloud use TLS, whichever transport is active (WiFi, Ethernet or 5G).

2. Prerequisites

- An active Corsano HCP account (cloud email + password).
- A study set up in the cloud for your site — the hub is registered into it.
- A Windows 10/11 PC with Internet access and a free USB port. No driver needed.
- The hub and a data-capable USB-C cable (not a charge-only cable).
- A network for the hub: 2.4 GHz WiFi credentials (5 GHz is not supported yet), an Ethernet outlet, or the built-in 5G SIM.
- 287-2B bracelets on firmware 8.36 or later.

Tip. Most failed onboardings come down to a charge-only USB cable, a 5 GHz-only WiFi network, or no study selected before clicking Register hub & provision.



3. Onboarding step by step

3.1 Step 1 — Sign in

Launch the console and sign in with your Corsano cloud credentials (the same as the HCP portal).

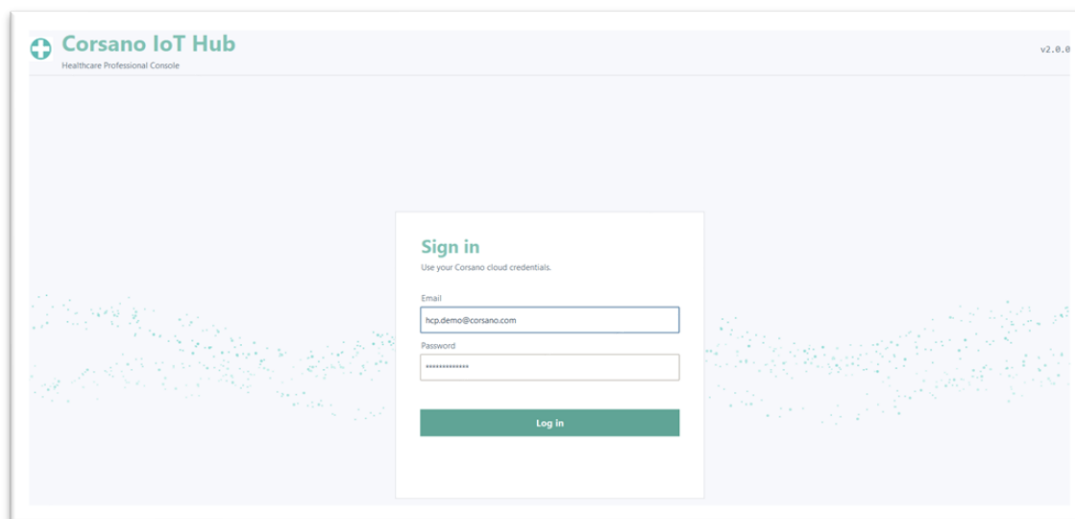


Figure 2 — Sign-in screen.

3.2 Step 2 — The main screen

Your studies are listed on the left. The right pane is the gateway panel — it reads “Plug a Corsano gateway via USB.” until a hub is connected.

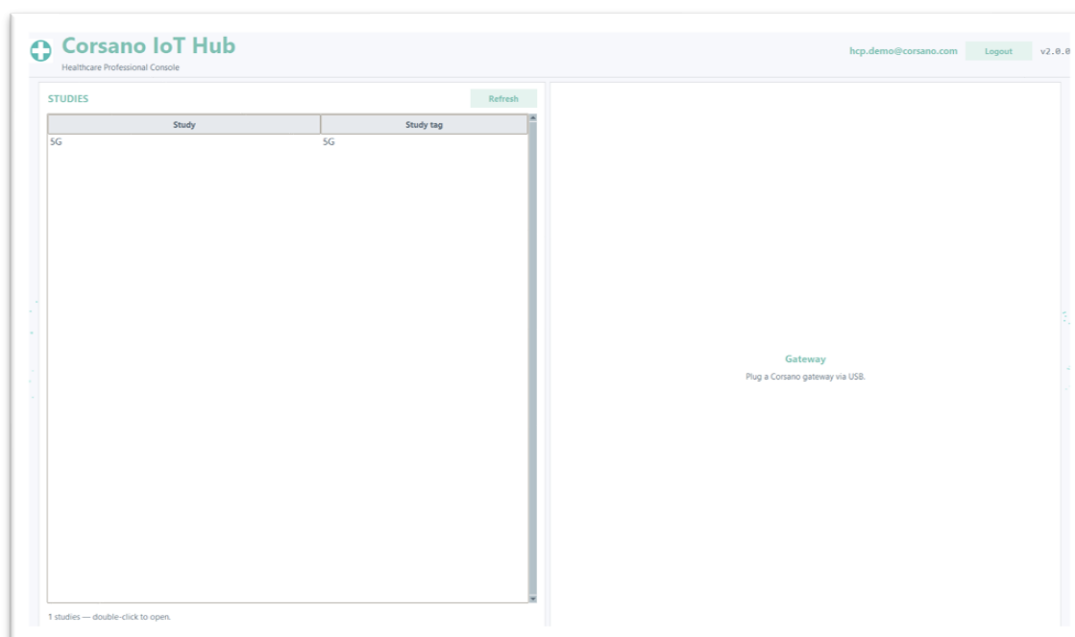


Figure 3 — Main screen, no hub connected yet.



3.3 Step 3 — Plug in the hub

Connect the hub over USB-C. Detection is automatic: within a second the GATEWAY banner shows its serial number and firmware. If the console offers a firmware update, apply it before continuing.

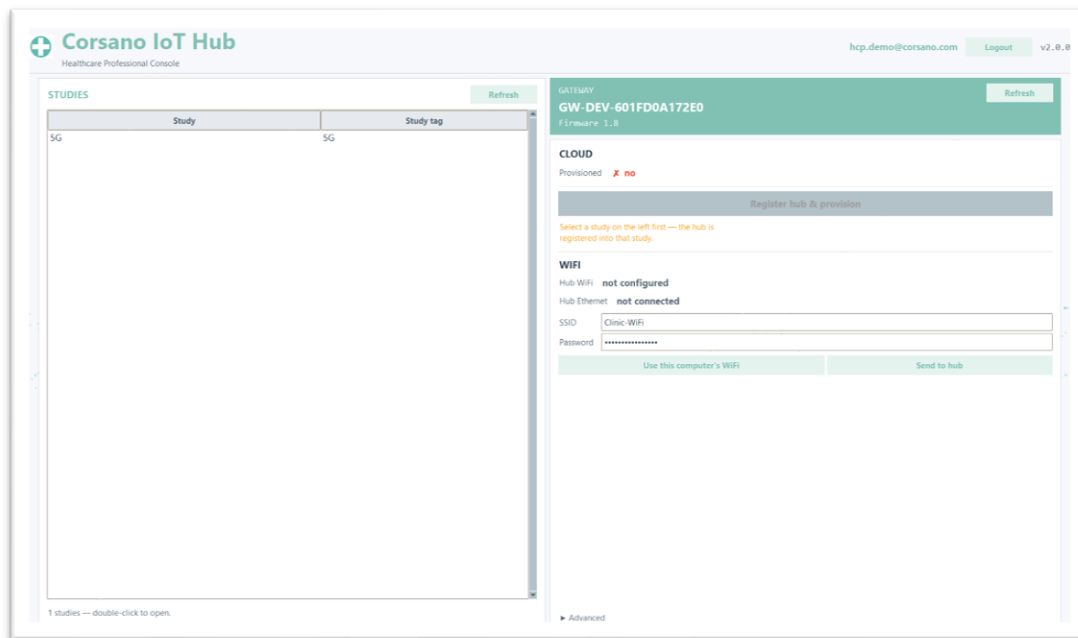


Figure 4 — Hub detected. The provision button stays greyed out until a study is selected.

3.4 Step 4 — Select the study

Double-click the target study. The left pane shows its sampling parameters (read-only) and Register hub & provision becomes active. A hub that was registered before gets its study pre-selected automatically.

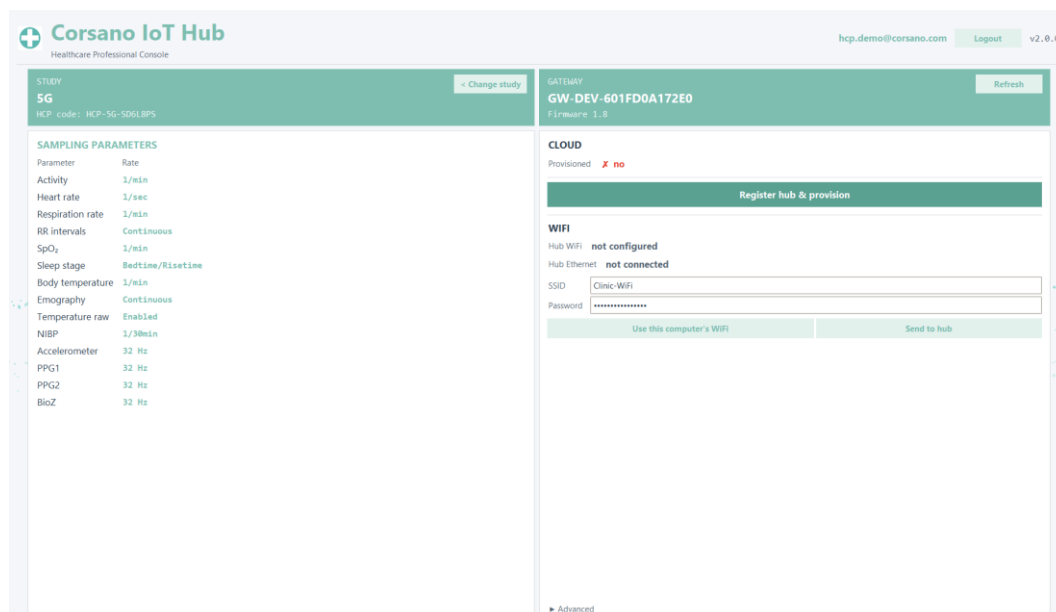


Figure 5 — Study selected — “Register hub & provision” is now active.



3.5 Step 5 — Enter the site WiFi (recommended)

In the WIFI section, enter the SSID and password of the network the hub should use — or click “Use this computer's WiFi” to copy the PC's current network. The credentials are pushed to the hub during provisioning. A hub on Ethernet or 5G can be provisioned with the fields left empty.

3.6 Step 6 — Register hub & provision

Click Register hub & provision. The console registers the hub in the Corsano cloud under your account, into the selected study, and configures the hub over USB. After a few seconds:

“GW-DEV-601FD0A172E0 is registered and provisioned. The hub is now fetching its bracelet list from the cloud. Assign bracelets to this gateway in the Corsano dashboard.”

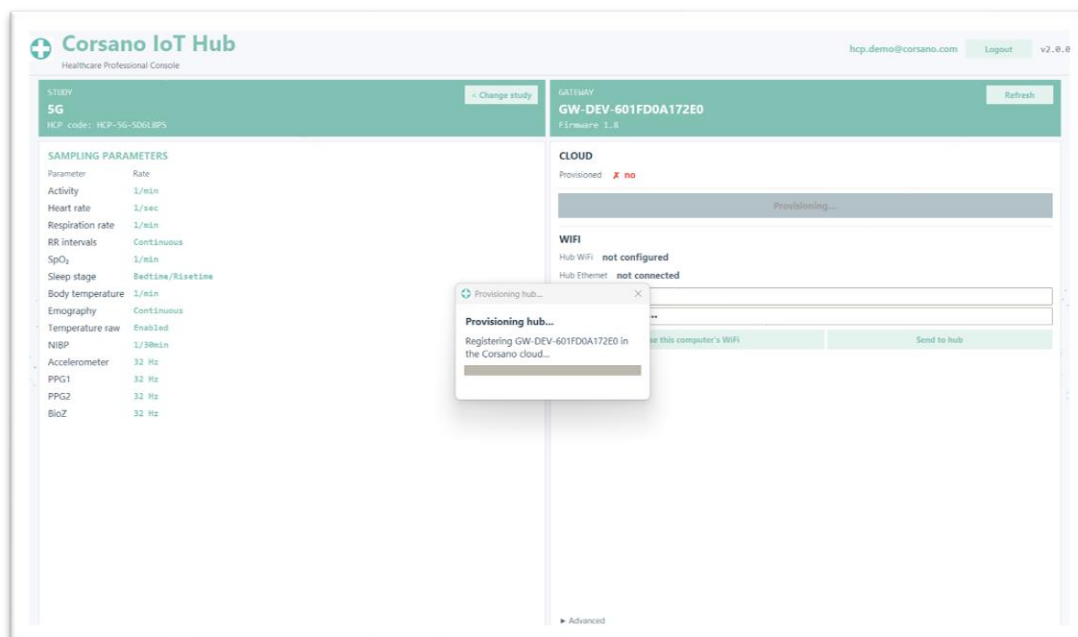


Figure 6 — Provisioning in progress.

3.7 Step 7 — Check the result

The CLOUD section now shows Provisioned ✓ yes and Cloud reachable ✓ yes — the hub reached the cloud over its own network (the WIFI section shows which one, with its IP). “Bracelet list: 0 bracelet(s)” is normal at this point.

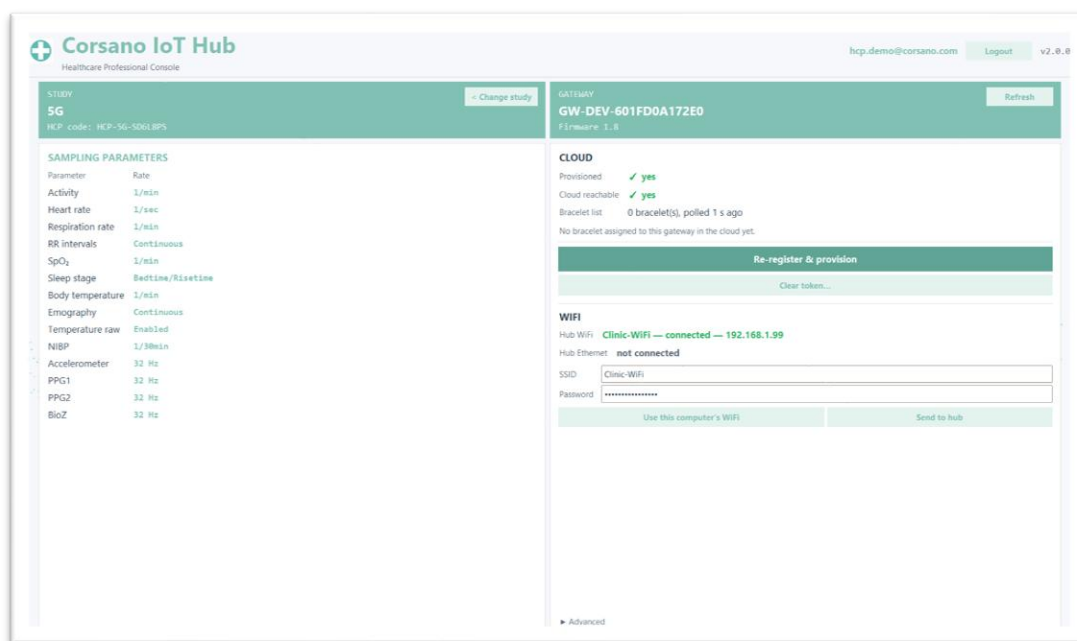


Figure 7 — Hub provisioned and online, no bracelets yet.

3.8 Step 8 — Assign the bracelets

Bracelets are assigned to the hub in the Corsano cloud dashboard, not in the console. The hub picks up the list within a minute — click Refresh to update the view.

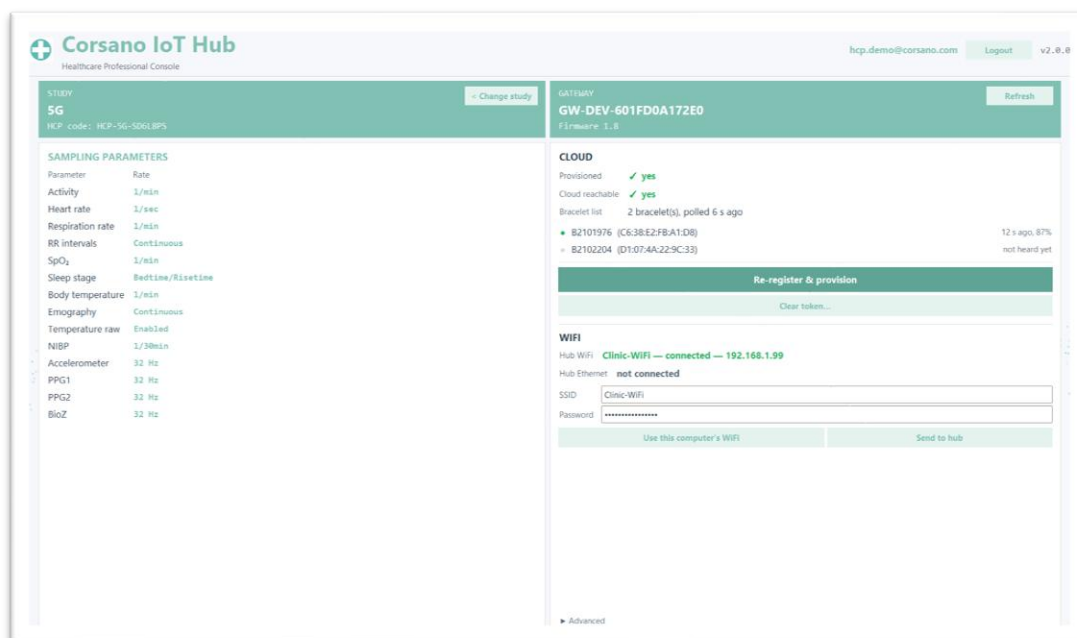


Figure 8 — Bracelet list filled in.



Display	Meaning
B2101976 (C6:38:...) — 12 s ago, 87%	Bracelet heard over BLE 12 s ago, battery 87 %.
B2102204 (D1:07:...) — not heard yet	Assigned in the cloud but not received yet — out of range, off, or firmware < 8.36.

3.9 Step 9 — Done

Unplug the USB cable — the hub now runs autonomously. Within a few minutes, vitals from worn bracelets appear in the Corsano portal.

Note. While the console is talking to the hub over USB, the hub pauses BLE collection; it resumes about 30 seconds after the last USB command. Data not flowing while the console is open on the hub is expected.

4. End-to-end verification

Run the checks in order — each failure points at the layer to investigate:

#	Check	Where	Expected
1	Hub detected over USB	Right pane	GATEWAY banner with SN + firmware.
2	Hub provisioned	CLOUD section	Provisioned ✓ yes.
3	Hub online	CLOUD section	Cloud reachable ✓ yes (over the hub's own network).
4	Hub network	WIFI section	“Hub WiFi ... — connected — <IP>” or “Hub Ethernet: connected”.
5	Bracelets assigned	CLOUD section	Bracelet list: N bracelet(s), N > 0.
6	Bracelets heard	Bracelet list	Recent “... ago” + battery on each worn bracelet.
7	Data in the cloud	Corsano portal	Patient vitals updating.

Note. Checks 1–5 pass but 6 fails → BLE side (bracelet off, out of range, old firmware). Checks 1–2 pass but 3 fails → hub network side (WiFi, Ethernet, SIM).



5. Troubleshooting

Symptom	Likely cause	Action
“Plug a Corsano gateway via USB.” never goes away	Charge-only cable, or the COM port is held by another app	Swap cable/USB port; close other serial tools.
Register hub & provision greyed out	No study selected	Double-click the study on the left (step 4).
Provisioning fails with “update the hub firmware”	Hub firmware too old for this app	Update the hub firmware, then retry.
“Could not obtain the gateway token ...”	Account lacks rights on the study, or no Internet on the PC	Check the HCP account and the PC connection; retry.
Provisioned ✓ but Cloud reachable ✗	Hub has no network (wrong WiFi, 5 GHz, no Ethernet/SIM)	Check the Hub WiFi line; resend 2.4 GHz credentials via Send to hub, or plug Ethernet.
“Hub WiFi: <ssid> — not connected”	Wrong password, or 5 GHz-only network	Fix the password; use a 2.4 GHz SSID.
“Bracelet list: 0 bracelet(s)” for good	Bracelets not assigned to this hub in the cloud	Assign them in the Corsano dashboard, wait ≤ 60 s, Refresh.
Bracelet stuck on “not heard yet”	Bracelet off, out of BLE range, or firmware < 8.36	Move/charge the bracelet; check its firmware; check BLE is enabled (Advanced section).
All green but no data in the portal	Console still connected over USB, or simple latency	Unplug USB, wait a few minutes, check again.

6. Day-to-day operations

- Refresh (GATEWAY banner) — re-reads the full hub state.
- Re-register & provision — re-runs registration with a fresh token; use after a study change or a factory reset. Idempotent.
- Clear token... — de-provisions the hub; it stops scanning and uploading. For decommissioning or moving a hub to another site.
- Send to hub (WIFI section) — pushes new WiFi credentials to an already-provisioned hub.
- Advanced (chevron at the bottom) — 5G self-test, LED, reboot, factory reset.



7. Support

Corsano Support — nicolas@corsano.com · <https://corsano.com/contact-us/>

When opening a ticket, include:

- the hub serial number (GATEWAY banner, e.g. GW-DEV-601FD0A172E0);
- the hub firmware and app version (top right, e.g. v2.0.0);
- the serial numbers of the affected bracelets;
- the window of impact (UTC) and a screenshot of the right pane.

End of document.